

New!



EASY

narrow support profile
batteries only on the outside
suitable for air transport
choice of three battery sizes



BASIC

square support profile
battery only inside
option to hang accessories



DUO / TRIO

square support profile
battery inside
additional battery outside for
double/triple range



Main advantages of HURT-e

One device for wheelchair user's mobility.
We are the manufacturers of this device ourselves.
Very light and quality construction.
Perfectly storage = small dimensions (in the car trunk, below the bed, behind the wardrobe ...).
Easy removal from the wheelchair thanks to quick release attachment.
Great power and range.
Option to connect a second battery for greater range.
Turn the wheelchair on the spot by turning the handlebars by 90°.
Choice of three levels of power.
For fixed and folding frame wheelchairs (folding retained).
Reverse gear.
Included: installation of attachment to the wheelchair, choice of basic colours, charger, casing, rear light for wheelchair.

KVADRU

move the handlebars toward you to
accelerate / push them away to brake
battery inside
option of additional external battery

Current price list and battery range at www.hurt-e.cz.

hurt-e.cz

Tel.: +420 777 327 910
E-mail: info@hurt-e.cz



HURT s. r. o.
Bambousek 664
281 26 Týnec nad Labem
Czech Republic

Parameters

Wheel size: 10 inches.

Motor: 48 V / 500 W.

Batteries manufactured in the Czech Republic.

Batteries: 48 V (AIR / BASIC / DUO / TRIO) 5.7 / 8.7 / 17.4 / 22.7 Ah.

Charging time (AIR / BASIC / DUO / TRIO) 2 / 6 / 12 / 17 hours.

Maximum range 15 / 30 / 60 / 80 km (depending on many variables).

Weight (EASY / BASIC / DUO / TRIO) 12 / 12 / 15 / 17 kg.

Maximum speed 25 km/h according to legislation.

Design speed up to 40 km/h. Adjustment at your own risk with warranty restrictions.

Choice of three power levels.

Maximum load 120 kg.

Dimensions 90 / 50 / 30 cm (height of device / width of handlebars / length including handle).

Disc brake with two independent brakes.

Brake lever lock.


Display showing speed, distance traveled, battery capacity.

Programmable user parameters.

Option for further individual adjustments.

Accessories

CROSS tire

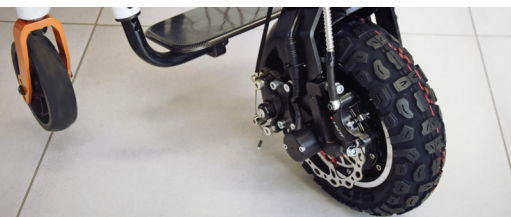

Easier passage through moderate terrain.

Multi-purpose bag



Easy installation on BASIC or KVADRU versions. We can send it by post.

Basket for handlebars



The metal basket can be easily removed thanks to the simple CLICK system. We can send by mail.

Hydraulic brake



Great brake effect in combination with a higher quality disc.

Left-hand accelerator



Left-handed throttle control.

Adjustable fork



For quadruplegic control. Adjustable.

Current price list and battery range at www.hurt-e.cz.